

INSERT WEAR GUIDE

Worn tips can diminish your scaling efficiency. Check your inserts **every week** to maintain optimal performance.

How the best perform

Hu-Friedy

#3 Beaver-tail

Original Propoxy Design

After Five® Right Low - Medium

After Five® Straight Low - Medium

After Five® Left Low - Medium

■ 25% Efficiency Loss (Reorder)

■ 50% Efficiency Loss (Discard)

#10 Universal Low - High

#100 Thin Low - Medium

#1000 Triple Bend Low - High

XT 25kHz | Low
30kHz | Low - Medium


XT Triple Bend Low - Medium

HOW TO USE YOUR WEAR GUIDE

1. Locate the diagram of the insert to be measured (use color as guide).
2. Position point **A** to point **B** as shown and/or align the curvature of the tip (circled) to match the guide.
3. Rotate insert until tip is flat against the card.
4. Evaluate tip condition, comparing length to red/blue bars.

Reorder tips when 25% efficiency loss (blue bar)

Discard tip when 50% efficiency loss (red bar)

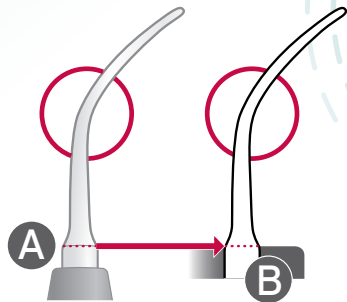

Learn more at [Hu-Friedy.com/Inserts](https://www.hu-friedy.com/inserts)

To reorder inserts call **1-800-Hu-Friedy** or your authorized Hu-Friedy dealer.

©2015 Hu-Friedy Mfg. Co., LLC. All rights reserved. HF-466/1015

Hu-Friedy