



## RUBBER DAM CLAMPS

REDUCE GLARE WITH SATIN STEEL®

How the best perform



# 48 WAYS TO INCREASE VISIBILITY



## FACT

Today's dental lighting is brighter and more intense than ever before. Yet most rubber dam clamps reduce visibility by reflecting this light back into your eyes, creating a blinding glare.

## CHALLENGE

Finding a rubber dam clamp that increases visibility and is available in a wide range of sizes and applications is often challenging.

## SOLUTION

Hu-Friedy's Satin Steel Rubber Dam Clamps offset bright clinical lighting. Hand-crafted to absorb and dissipate light instead of reflecting it back, these matte-finish rubber dam clamps are available in 48 sizes to help you perform at your best.

The result is clear, with Hu-Friedy's Satin Steel, you can experience a comfortable field of vision that's easier on the eyes.


SATIN STEEL RUBBER DAM CLAMP  
REFERENCE POSTER INSIDE.





# SATIN STEEL RUBBER DAM CLAMPS

*A complete line for anterior, premolar and molar applications*

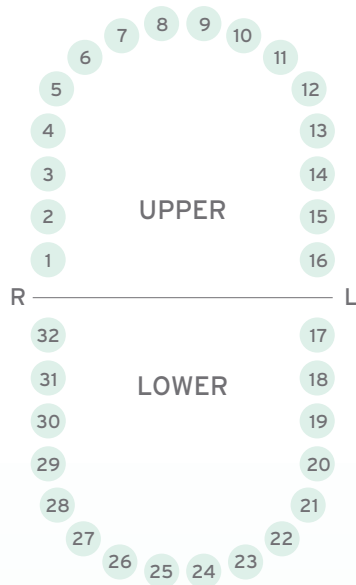

## ANTERIOR RUBBER DAM CLAMPS

**00**  
RDCM00

R L

Flat jaws and high bow for small lower canines and first premolars.

**6**  
RDCM6

R L

Slightly offset jaws for labial caries on anterior teeth.

**9**  
RDCM9

R L

Offset flat jaws and stiff bow for anterior teeth.

**9S**  
RDCM9S

R L

Offset flat serrated jaws for anterior teeth.

**212**  
RDCM212

R L

Flat jaws for labial caries on anterior teeth.


**212SA**  
RDCM212SA

R L

Slightly offset staggered jaws for odd shaped anteriors.

## ACCESSORIES

Hu-Friedy also has a full range of accessories useful for rubber dam clamp applications including forceps, punches, frames, templates, and an organizing board (shown on right).

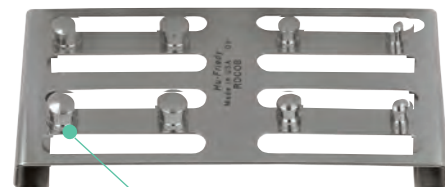

EIGHT PEGS SECURELY HOLD VARIOUS SIZED CLAMPS

How the best perform



## PREMOLAR RUBBER DAM CLAMPS

**0**  
RDCM0

Flat jaws and high bow for longer premolars.

This card shows the RDCM0 clamp, which has flat jaws and a high bow. It includes a top-down view of the clamp, a side view, and a diagram of the dental arch with green dots representing teeth. The diagram is labeled 'R' for right and 'L' for left.

**1**  
RDCM1

Slightly festooned jaws for upper premolars.

This card shows the RDCM1 clamp, which has slightly festooned jaws. It includes a top-down view of the clamp, a side view, and a diagram of the dental arch with green dots representing teeth. The diagram is labeled 'R' for right and 'L' for left.



SATIN STEEL PROVIDES SUPERIOR STRENGTH & CORROSION RESISTANCE

**1A**  
RDCM1A

Deeply festooned jaws for canines and premolars.

This card shows the RDCM1A clamp, which has deeply festooned jaws. It includes a top-down view of the clamp, a side view, and a diagram of the dental arch with green dots representing teeth. The diagram is labeled 'R' for right and 'L' for left.

**2**  
RDCM2

Flat jaws for premolars.

This card shows the RDCM2 clamp, which has flat jaws. It includes a top-down view of the clamp, a side view, and a diagram of the dental arch with green dots representing teeth. The diagram is labeled 'R' for right and 'L' for left.

**2A**  
RDCM2A

Flat jaws for larger premolars.

This card shows the RDCM2A clamp, which has flat jaws. It includes a top-down view of the clamp, a side view, and a diagram of the dental arch with green dots representing teeth. The diagram is labeled 'R' for right and 'L' for left.

**2AS**  
RDCM2AS

Flat serrated jaws for larger premolars.

This card shows the RDCM2AS clamp, which has flat serrated jaws. It includes a top-down view of the clamp, a side view, and a diagram of the dental arch with green dots representing teeth. The diagram is labeled 'R' for right and 'L' for left.

**W2**  
RDCMW2

Flat jaws for premolars.

This card shows the RDCMW2 clamp, which has flat jaws. It includes a top-down view of the clamp, a side view, and a diagram of the dental arch with green dots representing teeth. The diagram is labeled 'R' for right and 'L' for left.

**27N**  
RDCM27N

Slightly festooned jaws for premolars.

This card shows the RDCM27N clamp, which has slightly festooned jaws. It includes a top-down view of the clamp, a side view, and a diagram of the dental arch with green dots representing teeth. The diagram is labeled 'R' for right and 'L' for left.

**29**  
RDCM29

Flat jaws for lower premolars.


This card shows the RDCM29 clamp, which has flat jaws. It includes a top-down view of the clamp, a side view, and a diagram of the dental arch with green dots representing teeth. The diagram is labeled 'R' for right and 'L' for left.

**209**  
RDCM209

Small, flat jaws with stiff bow for lower premolars.

This card shows the RDCM209 clamp, which has small, flat jaws and a stiff bow. It includes a top-down view of the clamp, a side view, and a diagram of the dental arch with green dots representing teeth. The diagram is labeled 'R' for right and 'L' for left.

# MOLAR RUBBER DAM CLAMPS (GENERAL)

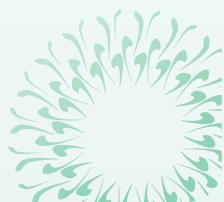

**8A**  
RDCM8A

Deeply festooned four-point jaws for molars.

**8AD**  
RDCM8AD

Deeply festooned jaws and extended bow for distal access to molars.

**W8A**  
RDCMW8A

Slightly festooned jaws for molars.

**10**  
RDCM10

Flat jaws designed for small upper left and lower right molars; paired with RDCM11.

**11**  
RDCM11

Flat jaws designed for small upper right and lower left molars; paired with RDCM10.

**14**  
RDCM14

Deeply festooned jaws for partially erupted or irregularly shaped molars.

**14A**  
RDCM14A

Deeply festooned jaws for partially erupted or irregularly shaped molars; larger than RDCM14.

**W14A**  
RDCMW14A

Deeply festooned jaws for partially erupted or irregularly shaped molars.

**18**  
RDCM18

Flat jaws and stiff bow for molars.

**24**  
RDCM24

Flat jaws for upper right and lower left molars; paired with RDCM25.

**25**  
RDCM25

Flat jaws for upper left and lower right molars; paired with RDCM24.

**27**  
RDCM27

Flat jaws for large molars.









**56**  
RDCM56

Slightly festooned jaws for large molars.











**56S**  
RDCM56S

Slightly festooned serrated jaws for large partially erupted or irregularly shaped molars.

## UPPER MOLAR RUBBER DAM CLAMPS

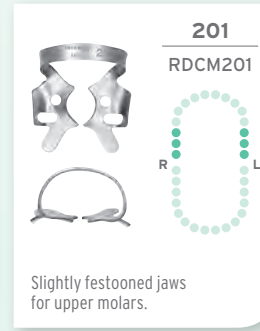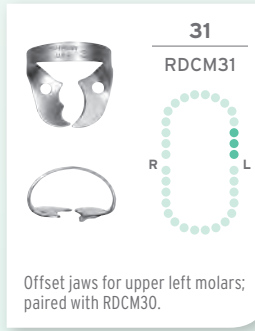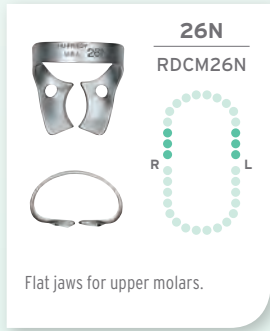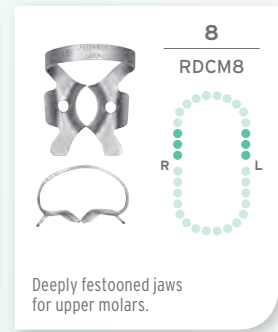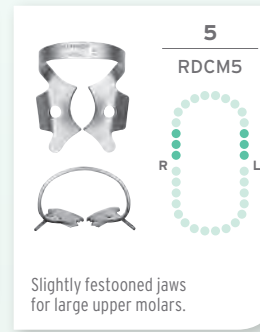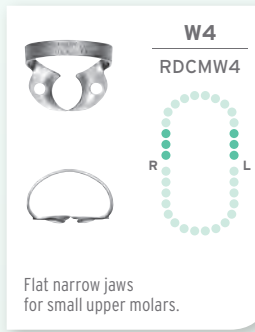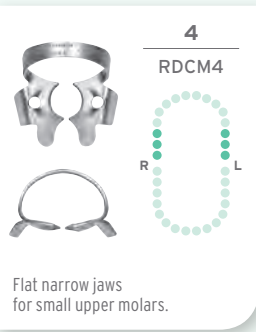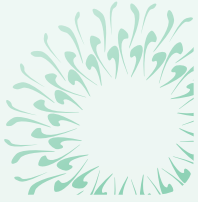

## LOWER MOLAR RUBBER DAM CLAMPS

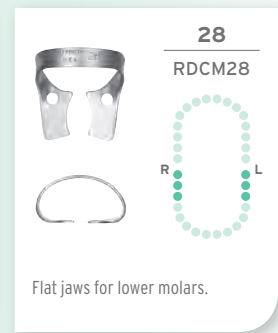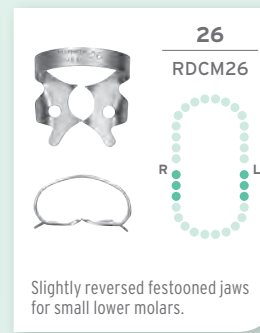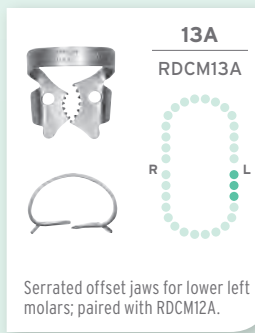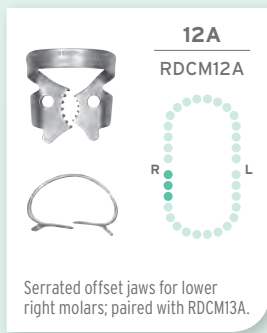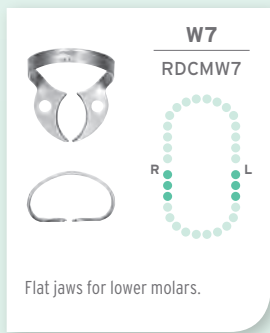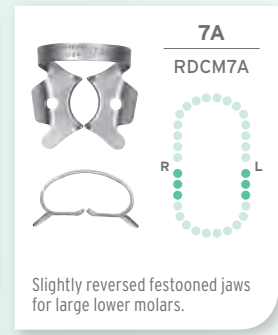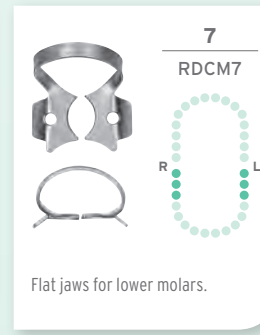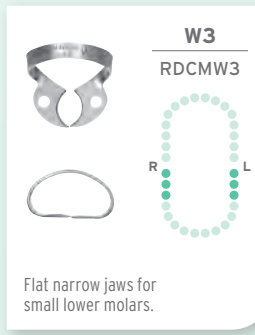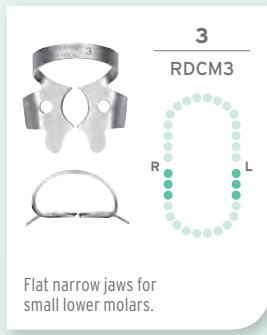
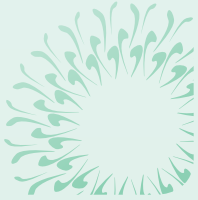


## RUBBER DAM KIT (RDKIT)

Basic items needed for rubber dam procedures conveniently packaged in a kit (also available individually).

DESCRIPTION	PART CODE
Rubber Dam Clamp Forceps	RDF
Rubber Dam Punch	RDP
Rubber Dam Frame-Adult	RDAF6
Universal Rubber Dam Template	RDTUNIV
FDI Rubber Dam Template	RDTFDI
Rubber Dam Clamp Organizing Board	RDCOB
#0 Satin Steel Rubber Dam Clamp	RDCM0
#1 Satin Steel Rubber Dam Clamp	RDCM1
#2 Satin Steel Rubber Dam Clamp	RDCM2
#7 Satin Steel Rubber Dam Clamp	RDCM7
#8 Satin Steel Rubber Dam Clamp	RDCM8
#8A Satin Steel Rubber Dam Clamp	RDCM8A
#14A Satin Steel Rubber Dam Clamp	RDCM14A
#212 Satin Steel Rubber Dam Clamp	RDCM212

## KITS AND SETS

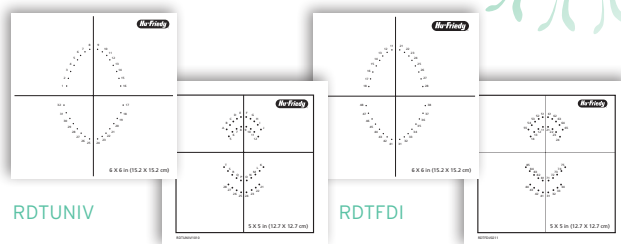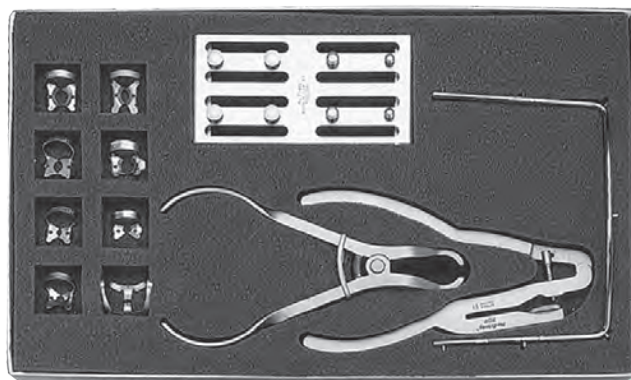

## RUBBER DAM CLAMP SET (RDCSET7)

A popular set for basic coverage. Other clamp sets also available.

DESCRIPTION	PART CODE
#00 Satin Steel Rubber Dam Clamp	RDCM00
#2 Satin Steel Rubber Dam Clamp	RDCM2
#7 Satin Steel Rubber Dam Clamp	RDCM7
#8 Satin Steel Rubber Dam Clamp	RDCM8
#8A Satin Steel Rubber Dam Clamp	RDCM8A
#9 Satin Steel Rubber Dam Clamp	RDCM9
#14A Satin Steel Rubber Dam Clamp	RDCM14A

# CROSS REFERENCE CHART

HU-FRIEDY PART CODE	COMPETITIVE CLAMP NUMBER	HU-FRIEDY PART CODE	COMPETITIVE CLAMP NUMBER	HU-FRIEDY PART CODE	COMPETITIVE CLAMP NUMBER
RDCM0	0	RDCMW7	W7	RDCM25	
RDCM00	00	RDCM8	8	RDCM26	26
RDCM1	1	RDCM8A	8A	RDCM26N	26N
RDCM1A	1A	RDCMW8A	W8A	RDCM27	27
RDCM2	2	RDCM8AD		RDCM27N	27N
RDCMW2	W2	RDCM9	9	RDCM28	28N
RDCM2A	2A	RDCM9S	9T	RDCM29	
RDCM2AS	2AT	RDCM10	10	RDCM30	
RDCM3	3	RDCM11	11	RDCM31	
RDCMW3	W3	RDCM12A	12A	RDCM56	56
RDCM4	4	RDCM13A	13A	RDCM56S	56T
RDCMW4	W4	RDCM14	14	RDCM201	
RDCM5	5	RDCM14A	14A	RDCM205	
RDCM6	6	RDCMW14A	W14A	RDCM209	209
RDCM7	7	RDCM18		RDCM212	212
RDCM7A	7A	RDCM24		RDCM212SA	212SA



For more information about  
Hu-Friedy's Total Solutions visit

**HU-FRIEDY.COM**


PERIODONTAL



RESTORATIVE



ORTHODONTIC

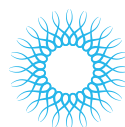


INSTRUMENT  
MANAGEMENT



ENDODONTIC



SURGICAL

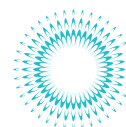

ULTRASONIC  
SCALING

