



ORTHO LINGUAL COLLECTION


REDEFINING THE ART
OF ORTHODONTICS

How the best perform




Hu-Friedy's

ORTHO LINGUAL COLLECTION

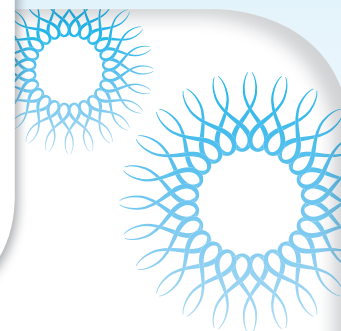

14% TO 29% THINNER
ORBIT FORMED JOINT
THAN STANDARD

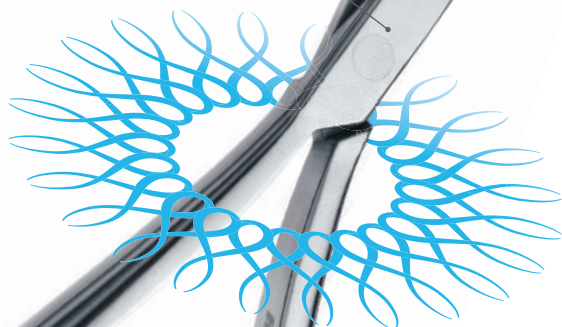

1.2" TO 2" LONGER
THAN STANDARD
HANDLE & .75" TO
1.2" LONGER THAN
LONG HANDLE

The Ortho Lingual Collection consists of ten instruments that feature:

- Optimal **intraoral access**
- Ideal **visibility** for clinician
- Enhanced **comfort** for patient
- More **time savings** allowing for additional patient education

Every ortho instrument from Hu-Friedy is created using the highest quality materials, processes and quality assurance.

DID YOU KNOW?

Hu-Friedy Orthodontic instruments are designed to withstand the rigors of steam sterilization while maintaining their beauty and functionality.

LINGUAL CUTTERS

Lingual Distal End Cutter

| 678-700

Wire size: .012" to .018"



Occlusal view of the Lingual Distal End Cutter cutting the excess archwire distal to the buccal tube.

Lingual Pin & Lig Cutter

| 678-701

Wire size: < = .012"




Occlusal view of the Lingual Pin & Lig Cutter cutting the excess steel ligature.

LINGUAL PLIERS

Lingual Weingart Pliers

| 678-702

Wire size: $\leq .030''$

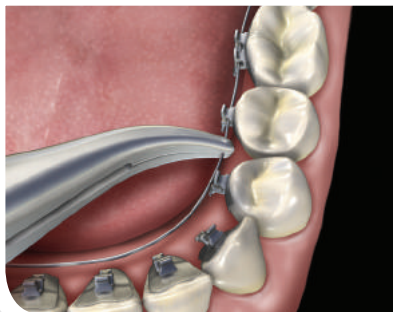

Occlusal view of the Lingual Weingart Pliers guiding the archwire into the bracket slot.

Lingual Weingart Pliers, 90° Bend

| 678-702-90

Wire size: $\leq .030''$



LINGUAL PLIERS



Lingual Mathieu Pliers

| 678-705




Occlusal view of the Lingual Mathieu Pliers tying steel ligature to hold the arch wire into bracket slot.



Lingual Bracket Removing Pliers

| 678-703

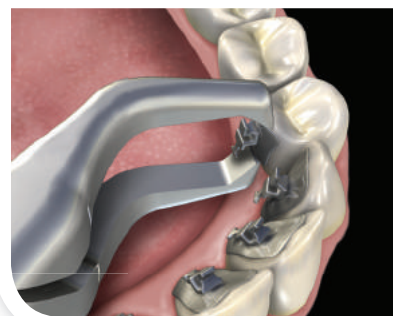


Lingual view of the Lingual Bracket Removing Pliers debonding the bracket.

LINGUAL PLIERS

Lingual Hammerhead® NiTi Pliers

| 678-704

Wire size: $\leq .025''$

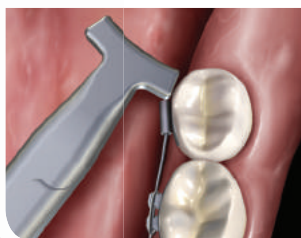

Occlusal view of the Lingual Hammerhead NiTi Pliers bending NiTi wire distal to the buccal tube.



Lingual Hammerhead NiTi Pliers with Serrations

| 678-704SER

Wire size: $\leq .025''$




Lingual Debonding Pliers

| 678-709



The arm hooks under the tie wing of the bracket while the other arm wedges under the bracket's adhesive interface. When squeezed, the wedged arm slides down the lingual surface of the tooth, successfully removing the bracket.


LINGUAL HAND INSTRUMENTS

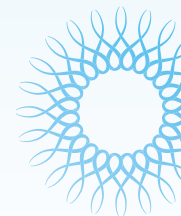

Lingual Ligature Director

| 678-706

Wire size: $\leq .021'' \times .025''$

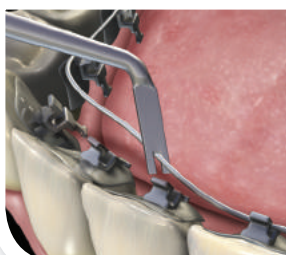

Ligature Side: Occlusal view of the Lingual Ligature Director guiding the archwire into the bracket slot.




Picker Side: Occlusal view of the Lingual Ligature Director tucking the excess steel ligature underneath the bracket.

Lingual Ligature Director, U Shape

| 678-707

Wire size: $\leq .021'' \times .025''$

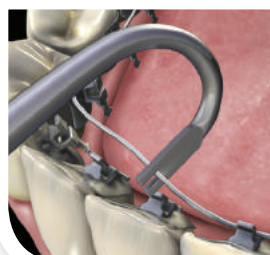


Occlusal view of the U Shaped Lingual Ligature Director guiding the archwire into the bracket slot.


Lingual Torquing Key with Picker

| 678-708

Slot size .018''

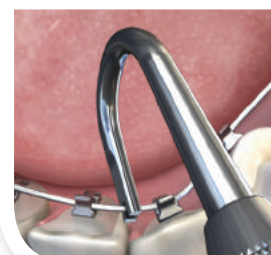

Occlusal view of the Torquing Key guiding the archwire into the bracket slot.

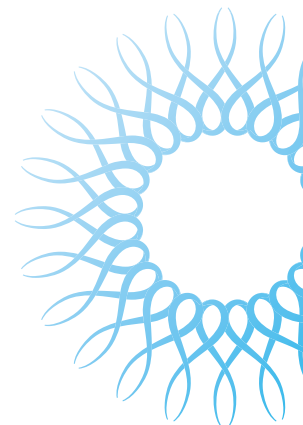

RECYCLE BROKEN & OLD INSTRUMENTS

HU-FRIEDY
PROGRAMS



When your instrument blades cannot be sharpened anymore, take advantage of Hu-Friedy's Envirodent® Instrument Recycling Program. Your old, broken instruments will be recycled to help create new car bumpers, bridges and street lights. For every 12 you recycle, we'll give you 1 FREE Hu-Friedy instrument up to \$47!



To start recycling visit
Envirodent.com

THE HU-FRIEDY ORTHO LINGUAL COLLECTION

CUTTERS

Lingual Distal End Cutter	678-700
Lingual Pin & Lig Cutter	678-701

PLIERS

Lingual Weingart Pliers	678-702
Lingual Weingart Pliers, 90°	678-702-90
Lingual Bracket Removing Pliers	678-703
Lingual Hammerhead NiTi Pliers	678-704
Lingual Hammerhead NiTi Pliers with Serrations	678-704SER
Lingual Mathieu Pliers	678-705
Lingual Debonding Pliers	678-709

HAND INSTRUMENTS

Lingual Ligature Director	678-706
Lingual Ligature Director, U Shape	678-707
Lingual Torquing Key with Picker	678-708

For more information about
Hu-Friedy's Total Solutions visit

HU-FRIEDY.COM


PERIODONTAL



RESTORATIVE



ORTHODONTIC

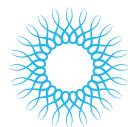


INSTRUMENT
MANAGEMENT



ENDODONTIC



SURGICAL

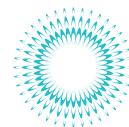

ULTRASONIC
SCALING



Hu-Friedy Mfg. Co., LLC: 3232 N. Rockwell St. | Chicago, IL 60618 | USA
Telephone: 1-800-Hu-Friedy (1-773-975-6100) | Website: Hu-Friedy.com

©2017 Hu-Friedy Mfg. Co., LLC. All rights reserved. [289] HF-899/1217

Hu-Friedy