

NEVI® SICKLE SCALERS

BETTER THAN EVER:
EVEREDGE 2.0

The Nevi Series features 4 scalers capable of scaling light to heavy calculus with minimal tissue trauma and superior adaption to allow precise access in difficult areas.

The Nevi instruments were specially designed with input from clinicians – dental hygienists like you who know: the details of oral anatomy, how the right bend or curve of the working end will allow you to successfully scale particular tooth surfaces, and that using the best instruments is what their patients deserve.


Each Nevi instrument was created from a combination of customer feedback and Hu-Friedy's tradition of quality, innovation and skill. The unique angles and distinct characteristics of each Nevi allow for access and reach into difficult areas, while our **EverEdge® 2.0 Technology** keeps the instruments sharper longer – adding up to a more efficient, comfortable scaling experience.



4 NEVI® POSTERIOR SCALER
SCNEVI49E2



NEVI® SICKLE SCALERS *Featuring EverEdge® 2.0 Technology*

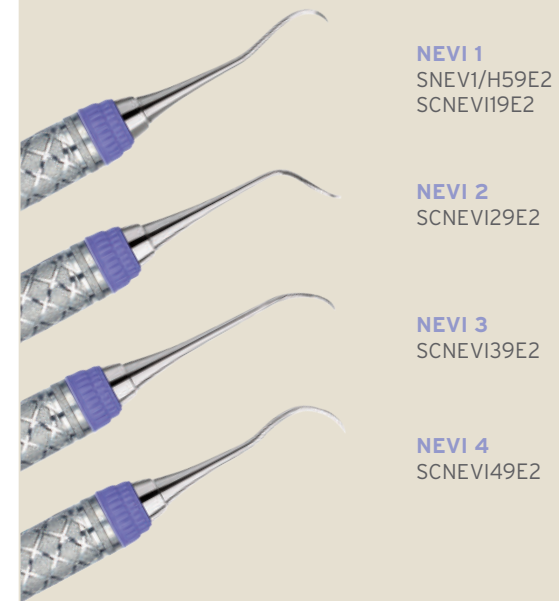

NEVI 1
SNEV1/H59E2
SCNEVI19E2

NEVI 2
SCNEVI29E2


NEVI 3
SCNEVI39E2

NEVI 4
SCNEVI49E2

NEVI AVAILABLE IN **RESIN 8 COLORS**

NEVI 1
SCNEVI1C8E2
NEVI 1/H5
SNEV1/H5C8E2
NEVI 2
SCNEVI2C8E2
NEVI 3
SCNEVI3C8E2
NEVI 4
SCNEVI4C8E2

For more information about
Hu-Friedy's Total Solutions visit
HU-FRIEDY.COM

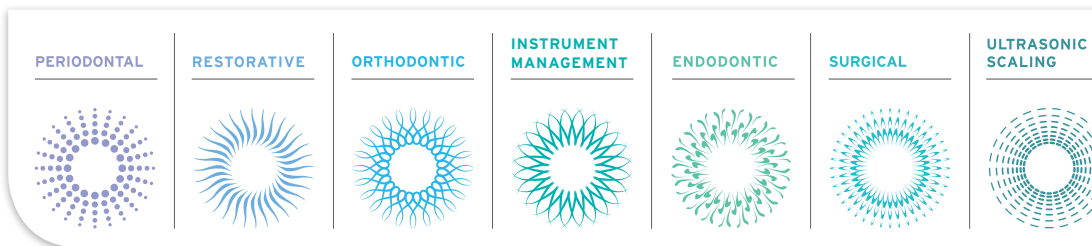

Hu-Friedy Mfg. Co., LLC: 3232 N. Rockwell St. | Chicago, IL 60618 | USA
Telephone: 1-800-Hu-Friedy (1-773-975-6100) | **Website:** Hu-Friedy.com

©2016 Hu-Friedy Mfg. Co., LLC. All rights reserved. HF-270/1116



NEVI® SICKLE SCALER COLLECTION

INNOVATIVE DESIGNS FOR
EFFECTIVE DEPOSIT REMOVAL



How the best perform



PERIODONTAL



ANTERIOR SICKLE SCALERS



NEVI® 1 | SCNEVI19E2



NEVI® 1 | SNEV1/H59E2



- Two distinct working ends provide solutions for anterior scaling needs
 - › **Sickle End:** A moderate length, thin scaler with two cutting edges and a blade that tapers to a point. Ideal for removing anterior interproximal calculus and stain.
 - › **Disc End:** A moderate sized, sharp rounded disc. Excellent for removal of anterior lingual stain.
- The Nevi 1 scaler is also available with a longer sickle working end



MESIAL LEFT
CENTRAL INCISOR



LINGUAL LEFT
LATERAL INCISOR



POSTERIOR SICKLE SCALERS

NEVI® 2 | SCNEVI29E2



- Perfect for **light to moderate** interproximal calculus removal
- The thinnest dual tapered blades in the Nevi collection remove calculus with minimal tissue trauma
- Deep contra-angles are paired with extra blade length for easy access under contacts



MESIAL RIGHT
FIRST MOLAR



DISTAL LEFT
SECOND MOLAR

NEVI® 3 | SCNEVI39E2



- Used for light to moderate calculus removal in **tight, narrow posterior interproximal areas**
- Ultra slender, dual cutting blades for effective placement
- Contra-angle bends in shank offer superior reach and adaptation to posterior interproximal areas



MESIAL LEFT
FIRST MOLAR



DISTAL LEFT
SECOND PREMOLAR

POSTERIOR SICKLE SCALERS

NEVI® 4 | SCNEVI49E2



- Ideal for scaling **heavy posterior interproximal deposits**
- Increased shank size and rigidity generates effective removal of tenacious, heavy and burnished deposits
- Working ends with dual tapered cutting edges allow for precise interproximal access
- Largest blade in the Nevi collection



MESIAL RIGHT
FIRST PREMOLAR



DISTAL RIGHT
SECOND MOLAR

