



AIR-FLOW® POWDERS

AIR-FLOW powders are designed for precise compatibility with AIR-FLOW devices in order to deliver the most comfortable and efficient experience for the patient and clinician. Whether your practice seeks cosmetic stain removal or a cutting edge solution for biofilm management and perio maintenance, Hu-Friedy EMS has the right powder for you.




NEW



ERYTHRITOL AIR-FLOW PLUS POWDER

Erythritol is the newest innovation in dental powders with a micron size of only 14. This ultrafine powder is a sugar alcohol that is already used in dentistry and the food industry. It is safe to use subgingivally and on/around delicate surfaces such as implants and restorations, but the harder material is twice as efficient as glycine at removing stain.

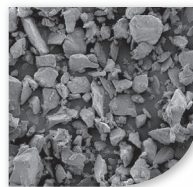

PARTICLE
Size



USAGE
Area*




LIGHT-MODERATE
Stain Removal



GLYCINE AIR-FLOW PERIO POWDER

Glycine is a low abrasive amino acid-based powder that is safe on almost all surfaces. This fine-grain powder is ideal for subgingival use and aiding in disease prevention for periodontitis, mucositis and peri-implantitis. AIR-FLOW Perio powder is 25 microns for the most gentle patient experience.

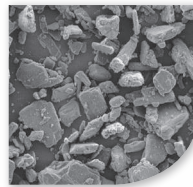

PARTICLE
Size



USAGE
Area*




LIGHT-MODERATE
Stain Removal



SODIUM BICARBONATE AIR-FLOW CLASSIC POWDERS *featuring NEW COMFORT FORMULA*

Sodium Bicarbonate is the original material used in air polishing powders. It is intended for heavy stain removal due to the grain structure and size. Most sodium bicarbonate powders range from 60-120 microns, but AIR-FLOW Classic's new formula is just 40 microns for equivalent efficiency and efficacy with added comfort.

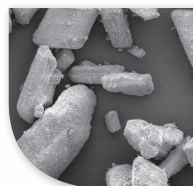

PARTICLE
Size




USAGE
Area



HEAVY
Stain Removal



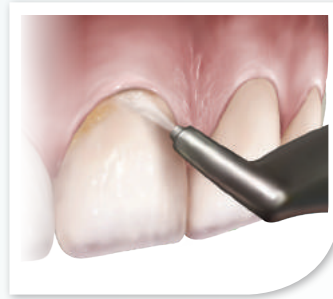

WHICH POWDER IS RIGHT FOR EACH PROCEDURE?

	CLASSIC	PERIO	PLUS
Heavy Stain	•		
Light-Moderate Stain	•	•	•
Subgingival Biofilm		•	•
Implants		•	•
Ortho		•	•
Surface Prep		•	•
Sealant Replacement		•	•
Dentures		•	•
Perio Maintenance		•	•
Cavitated Areas		•	•
Restorative Materials		•	•

“ I use Plus powder for minimally-invasive biofilm removal on enamel, exposed dentin, inside deep periodontal pockets, on implants, around orthodontic bands, and on delicate restorative materials – *everywhere!* It is also effective on stain removal, saving me substantial time while protecting my patients’ teeth. No wonder my patients prefer this gentle treatment compared to traditional methods! ”

–KAREN DAVIS, RDH, BSDH



AVAILABLE POWDERS		POWDER COMPATIBILITY					
IN PACKS OF 4	PART CODE	AIR-FLOW handy 3.0	AIR-FLOW handy 3.0 Perio/Premium	AIR-FLOW Master®	AIR-FLOW Master Piezon®	AIR-FLOW S1	AIR-FLOW S2
AIR-FLOW Classic Powder, Lemon	DV-1001	✓		✓	✓	✓	✓
AIR-FLOW Classic Powder, Mint	DV-1010	✓		✓	✓	✓	✓
AIR-FLOW Classic Powder, Cherry	DV-1006	✓		✓	✓	✓	✓
AIR-FLOW Classic Powder, Neutral	DV-1009	✓		✓	✓	✓	✓
AIR-FLOW Perio Powder	DV-110/A/HS	◆	✓	✓	✓		
AIR-FLOW Plus Powder	DV-117A	◆	✓	✓	✓		

✓ optimal ◆ acceptable

Learn more at Hu-Friedy.com/Biofilm

Hu-Friedy Mfg. Co., LLC: 3232 N. Rockwell St. | Chicago, IL 60618 | USA
Telephone: 1-800-Hu-Friedy or 1-773-975-6100 | Website: Hu-Friedy.com

©2018 Hu-Friedy Mfg. Co., LLC. All rights reserved. [557_2] HFL-087/0518

