

PERIODONTAL



## A SHARPER EDGE FOR BETTER CARE

| EverEdge® Technology

How the best perform



## EVEREDGE® TECHNOLOGY

Hu-Friedy's EverEdge Technology is unlike anything you've experienced before in a scaler. We've combined state-of-the-art technology in metallurgy, heat treatment and cryogenics to create a superior stainless steel that keeps your scalers and currettes sharper longer. That means less frequent sharpening, less hand fatigue, and greater comfort throughout the day.

“ As an educator, I always recommend Hu-Friedy's EverEdge instruments to my dental hygiene students. When well maintained, they can be sharpened so quickly and easily. A real advantage of the fine sharp cutting edge is **less hand fatigue for the hygienist and better comfort for the patient**, so both the clinician and patient can feel the benefits of EverEdge. ”

— Jill Nield Gehrig, RDH, MA  
Dean Emeritus, Division of Allied Health

### STAYS SHARPER LONGER

The revolutionary stainless steel alloy is super-durable and stays sharper longer so the clinician can spend less time sharpening.

### PROPRIETARY PROCESS

Our proprietary processing ensures that the superior edge retention and wear characteristics of EverEdge Technology will last the entire life of the instrument. Since it's not a superficial coating, they can be resharpened like any other scaler to extend the useful life.

### FINELY-HONED FINISH

EverEdge Scalers feature a hand-crafted finish from Hu-Friedy's professional artisans that keeps edges effectively sharp and allows for efficient scaling.

### DIAMOND KNURL


The unique diamond knurl pattern on the large diameter handle provides a confident and comfortable grasp. The #9 EverEdge handle was also designed to be the optimal weight for reduced hand fatigue.

### COMFORT ZONE CONE WITH SIGNATURE SERIES GRIP

The cushioned grip provides a smooth transition and increased comfort. With distinctive Ocean Blue grips, EverEdge scalers can be easily identified in a set up or customized for your office's needs.

### NOW AVAILABLE IN RESIN 8 COLORS

These lightweight resin handles use a color-coded system to easily identify scalers and currettes that correspond to anatomical areas of the mouth.

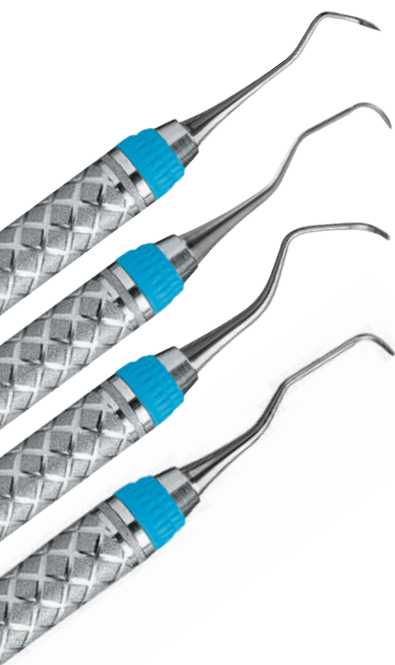

## SCALER CHECK-UP


Hu-Friedy uses the highest quality materials and processes to create EverEdge scalers' long lasting cutting edges but they still need to be regularly maintained and then replaced after a certain amount of time. The useful life of a scaler depends on several factors including: frequency of use, frequency of sharpening, sterilization method and clinical applications. Hu-Friedy recommends that you regularly audit your scalers to determine when they should be replaced.

### AUDIT CHECKLIST:

- Retipped instruments
- Dull/chipped cutting edge or point
- Broken or bent tips
- Thin/worn-out blades
- Pointed curettes
- Corrosion

### DID YOU KNOW?

Curettes are designed with a rounded toe for patient comfort during subgingival scaling. Sharpening can turn this "toe" into a point and cause patient discomfort when used subgingivally.

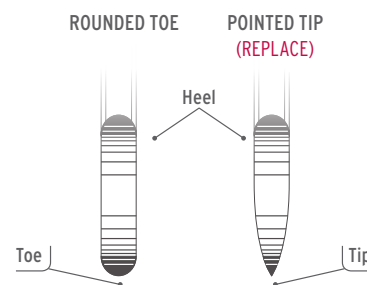

## TOP 10 SCALERS

1	H6/H7	SH6/79
2	204S	S204S9
3	Nevi 4	SCNEV149
4	Gracey 13/14	SG13/149
5	Gracey 11/12	SG11/129
6	Columbia 13/14	SC13/149
7	Barnhart 5/6	SBH5/69
8	Gracey 1/2	SG1/29
9	Nevi 2	SCNEV129
10	H5/33 Jacquette	SH5/339

For more information about Hu-Friedy's Total Solutions visit

**HU-FRIEDY.COM**


PERIODONTAL



RESTORATIVE



ORTHODONTIC

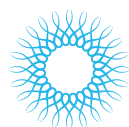


INSTRUMENT MANAGEMENT



ENDODONTIC



SURGICAL

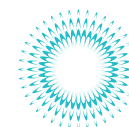

ULTRASONIC SCALING



Hu-Friedy Mfg. Co., LLC: 3232 N. Rockwell St. | Chicago, IL 60618 | USA  
Telephone: 1-800-Hu-Friedy (1-773-975-6100) | Website: Hu-Friedy.com

©2014 Hu-Friedy Mfg. Co., LLC. All rights reserved. HF-479/0914

**Hu-Friedy**